

NF κ B p65 Polyclonal Antibody

CATALOG #:	3038-100
AMOUNT:	100 μ g
LOT #:	_____
HOST:	Rabbit
IMMUNOGEN:	Synthetic peptide surrounding amino acid 274 of mouse NF κ B p65 (ID# BV-129)

FORMULATION:

100 μ g (0.2 mg/ml) immunoaffinity purified rabbit anti-NF κ B polyclonal antibody in phosphate buffered saline (PBS), pH 7.2, containing 25% glycerol, 1% BSA, 0.02% thimerosal.

SPECIES REACTIVITY: human, mouse, rat

STORAGE CONDITIONS:

Store at -20°C. For long-term storage, aliquot and freeze at -70°C. Avoid repeated freeze/thaw cycles.

BACKGROUND:


Nuclear factor kappa B (NF κ B) was identified as a sequence specific transcriptional activator that binds to the intronic enhancer of kappa light chain gene in B lymphocytes. NF κ B is a heterodimer that consists of a 50 kDa DNA binding subunit (p50) and a 65 kDa transactivation subunit (p65/RelA). Both of these subunits exhibit sequence homology to the protooncogene c-Rel. The p50 has an isoform called p49/p52, and both proteins are derived from the amino-terminal of precursor protein p105 and p100. The p50/p65 heterodimer remains in the cytosol in an inactive form as a complex with its inhibitor, I κ B. Upon stimulation of cells by a wide variety of stimuli such as lipopolysaccharide (LPS), pro-inflammatory cytokines (IL-1 & TNF, etc.), and viral infection, I κ B is phosphorylated and degraded by proteasome. The active NF κ B heterodimer is translocated into the nucleus and induces gene expression.

SPECIFICITY:

Detects NF κ B p65 and p50 (in less extent) in samples of human, mouse, and rat origins.

APPLICATIONS:

The antibody can be used for Western blotting (1-2 μ g/ml) and immunocytochemistry (10-20 μ g/ml). However, the optimal conditions should be determined individually.

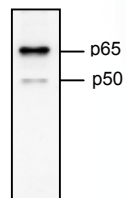

Western blot analysis of Jurkat cells using NF- κ B antibody.

RELATED PRODUCTS:

Apoptosis Detection Kits & Reagents

- Annexin V Kits & Bulk Reagents
- Caspase Assay Kits & Reagents
- Mitochondrial Apoptosis Kits & Reagents
- Nuclear Apoptosis Kits & Reagents
- Apoptosis Inducers and Set
- Apoptosis siRNA Vectors

Cell Fractionation System

- Mitochondria/Cytosol Fractionation Kit
- Nuclear/Cytosol Fractionation Kit
- Membrane Protein Extraction Kit
- Cytosol/Particulate Rapid Separation Kit
- Mammalian Cell Extraction Kit
- FractionPREP Fractionation System

Cell Proliferation & Senescence

- Quick Cell Proliferation Assay Kit
- Senescence Detection Kit
- High Throughput Apoptosis/Cell Viability Assay Kits
- LDH-Cytotoxicity Assay Kit
- Bioluminescence Cytotoxicity Assay Kit
- Live/Dead Cell Staining Kit

Cell Damage & Repair

- HDAC Fluorometric & Colorimetric Assays & Drug Discovery Kits
- HAT Colorimetric Assay Kit & Reagents
- DNA Damage Quantification Kit
- Glutathione & Nitric Oxide Fluorometric & Colorimetric Assay Kits

Signal Transduction

- cAMP & cGMP Assay Kits
- Akt & JNK Activity Assay Kits
- Beta-Secretase Activity Assay Kit

Adipocyte & Lipid Transfer

- Recombinant Adiponectin, Survivin, & Leptin
- CETP Activity Assay & Drug Discovery Kits
- PLTP Activity Assay & Drug Discovery Kits
- Total Cholesterol Quantification Kit

Molecular Biology & Reporter Assays

- siRNA Vectors
- Cloning Insert Quick Screening Kit
- Mitochondrial & Genomic DNA Isolation Kits
- 5 Minutes DNA Ligation Kit
- 20 Minutes Gel Staining/Destaining Kit
- β -Galactosidase Staining Kit & Luciferase Reporter Assay Kit

Growth Factors and Cytokines

Monoclonal and Polyclonal Antibodies