

KR104 Anti Mouse DC-STAMP Polyclonal Antibody

Package Size: 25µg (100µL / vial)

Format Rabbit polyclonal antibody: 0.25mg/mL

Buffer: PBS [containing 2% Block Ace as a stabilizer, 0.1% Proclin as a bacteriostat]


Storage Store: below -20°C

Once thawed, store at 4°C. Repeated freeze-thaw cycles should be avoided.


Purification method: This antibody was purified from rabbit serum immunized with synthesized partial peptide of mouse DC-STAMP(FIND) by peptide affinity chromatography.

Working dilution: For Immunohistochemistry ; 4µg/mL

① Anti DC-STAMP antibody

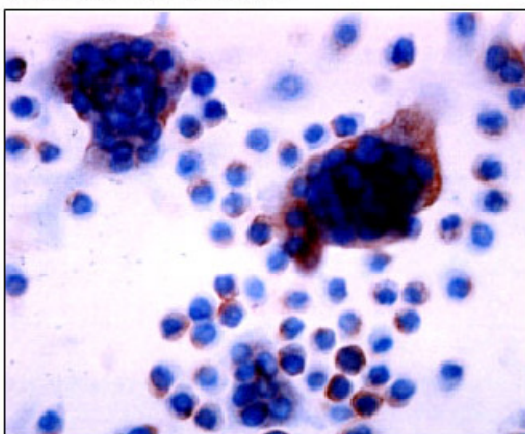

② Control rabbit IgG⁺

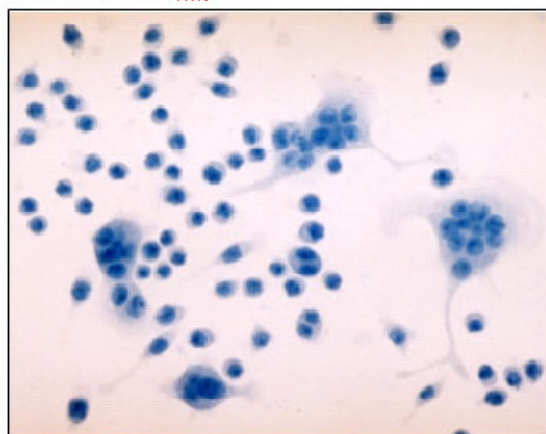

Immunohistochemistry

Sample :

RAW-D cells stimulated with RANKL+TNF-alpha

1. Anti DC-STAMP antibody
2. Control rabbit IgG

□Reference□

1. Hartgers F. C. et al. : Eur. J. Immunol. 2000 Dec;30(12):3585-90.
2. Hartgers F. C. et al. : Immunogenetics. 2001 Mar;53(2):145-9.
3. Staeger H. et al. : Immunogenetics. 2001 Mar;53(2):105-13.
4. Rho J. et al.: DNA Cell Biol. 2002 Aug;21(8):541-9.
5. Kukita T, et al. : J Exp Med. 2004 Oct 4;200(7):941-6. *
- 6 Yagi M. et al.: J Exp Med. 2005 Aug 1;202(3):345-51.

Gentaur Europe

Avenue de l'Armée 68
B-1040 Brussels
Tel: +32 16 58 90 45
Fax: +32 16 50 90 45
info@gentaur.com
www.gentaur.com

Gentaur Germany

IME, Forckenbeckstrasse 6
D-52074 Aachen Germany
Tel: 0241 6085 13140
Fax: 0241 6085 33033
info@bioprice.com
www.bioprice.com

Gentaur France

9, Rue Lagrange
75005 Paris, France
Tel: +33(0)143 25 0150
Fax: +33(0)143 25 0160
info@clonagen.com
www.clonagen.com

Gentaur Italy

Tel: 02 36 00 65 93
Fax: 02 36 00 65 94
info@labprice.com
www.labprice.com